

## Medical & Laboratory cold storage system Ultra-Low Temperature Freezer Temperature -60/ -86°C Capacity: 50~936L

### Handle

Handle for easy operation, door lock for safety use.



### MicrocomputerControl

Digital temperature display, ensure precise and stable running. USB port to read temperature and alarm conditions in recent 3 months. Able to be change to touch screen display.

### USB & Remote Alarm Port

USB port make it easy to read and store temperature record, remote alarm port to connect remote alarm system.



### Heavy Duty Castors

Castors with break, easy moving and stop.



### 3 Kinds of Optional Cooling System

1. Cascade system with refrigerant R404a and R507a, most widely used system, easy to repair.
2. Cascade system with refrigerant R290 and R170, green refrigerant.
3. One compressor with mixed refrigerant. Energy saving and better price.

Up to 43°C ambient temperature, high efficient freezing system and good thermal insulation properties

Low noise design  
designNoise≤50db

All 304 stainless steel inner door, with foaming layer and small handle.

4 layers door seal made of high quality silicone






ULT freezer is specially designed and manufactured for long term storage of various biological products, including viruses, germs, erythrocytes, leucocytes and cutis. Applications can be found in blood banks, hospitals, epidemic prevention services, and research institutes, laboratories in electronic and chemical plants, biological engineering institutes and marine fishery companies.

Temperature profile for -60°C



Temperature profile for -86°C

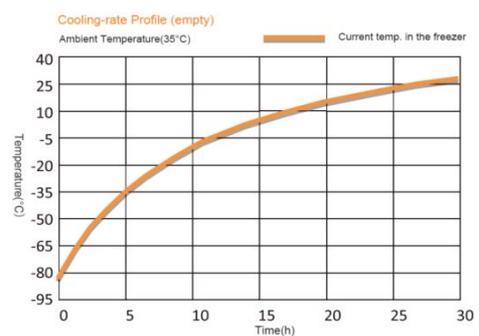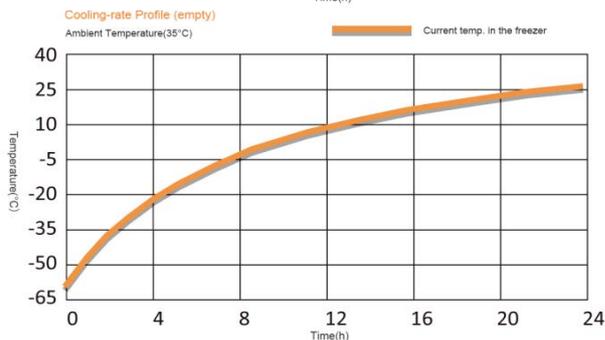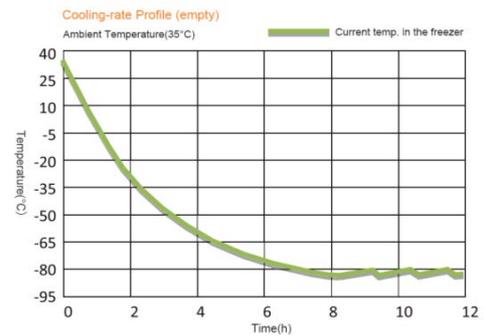

## ST-MDF-86V50



Model	ST-MDF-86V50	Capacity(L)	50
External Size(mm)	540 X 620 X 1115	N.W./G.W.(kg)	75/95
Inner Size(mm)	361 X 429 X 406		
Package Size(mm)	690 X 760 X 1170	Insulation	PURF
Voltage	220V; 50/60Hz	External Material	Color sprayed steel
Power(W)	450	Inner Material	304 Stainless steel
Setting Temp.(°C)	-40~-86	Shelves	1
Ambient Temp. (°C)	10-32	Noise(dBA)	≤49.5
Alarm System	High/low temperature, sensor failure		
Optional Accessories	Remote alarm system, cold storage racks & boxes		

## ST-MDF-86V188E



Model	ST-MDF-86V188E	Capacity(L)	188
External Size(mm)	830 X 810 X 1607.5	N.W./G.W.(kg)	150/195.5
Inner Size(mm)	465 X 465 X 830		
Package Size(mm)	880 X 905 X 1820	Insulation	PURF
Voltage	220V; 50/60Hz	External Material	Color sprayed steel
Power(W)	600	Inner Material	304 Stainless steel
Setting Temp.(°C)	-40~-86	Shelves	2
Ambient Temp. (°C)	10-32	Noise(dBA)	≤51
Alarm System	High/low temperature, sensor failure, controller		
Optional Accessories	Cold storage racks & boxes, chart recorder		

## ST-MDF-86V340E



Model	ST-MDF-86V340E	Capacity(L)	340
External Size(mm)	855 X 1006 X 1920	N.W./G.W.(kg)	181.5/234
Inner Size(mm)	488 X 600 X 1140		
Package Size(mm)	935 X 1070 X 2135	Insulation	PURF
Voltage	220V; 50/60Hz	External Material	Color sprayed steel
Power(W)	1350	Inner Material	304 Stainless steel
Setting Temp.(°C)	-40~-86	Shelves	3
Ambient Temp. (°C)	10-38	Noise(dBA)	≤50
Alarm System	High/low temperature, sensor failure, door ajar		
Optional Accessories	Remote alarm system, cold storage racks & boxes, chart recorder		

## ST-MDF-86V340



Model	ST-MDF-86V340	Capacity(L)	340
External Size(mm)	854 X 1072 X 1907	N.W./G.W.(kg)	249.5/302
Inner Size(mm)	488 X 607 X 1140		
Package Size(mm)	935 X 1070 X 2135	Insulation	PURF
Voltage	220V; 50/60Hz	External Material	Color sprayed steel
Power(W)	1350	Inner Material	304 Stainless steel
Setting Temp.(°C)	-60~-86	Shelves	3
Ambient Temp. (°C)	10-38	Noise(dBA)	≤50
Alarm System	High/low temperature, sensor failure, door ajar, power failure(72h), low battery		
Optional Accessories	Remote alarm system, cold storage racks & boxes, chart recorder		

## ST-MDF-86V588(VIP)



Model	ST-MDF-86V588	Capacity(L)	588
External Size(mm)	872 X 1192 X 1994.5	N.W./G.W.(kg)	249/310.5
Inner Size(mm)	606 X 738 X 1310		
Package Size(mm)	974 X 1118 X 2183	Insulation	PURF
Voltage	220V; 50/60Hz	External Material	Color sprayed steel
Power(W)	815	Inner Material	304 Stainless steel
Setting Temp.(°C)	-60~-86	Shelves	3
Ambient Temp. (°C)	10-32	Noise(dBA)	≤53
Alarm System	High/low temperature, sensor failure, door ajar, power failure(72h), low battery		
Optional Accessories	Remote alarm system, cold storage racks & boxes, chart recorder		

## ST-MDF-86V688



Model	ST-MDF-86V688	Capacity(L)	688
External Size(mm)	1354 X 1072 X 1907	N.W./G.W.(kg)	367/446
Inner Size(mm)	990 X 607 X 1140		
Package Size(mm)	1475 X 1090 X 2145	Insulation	PURF
Voltage	220V; 50/60Hz	External Material	Color sprayed steel
Power(W)	1540	Inner Material	304 Stainless steel
Setting Temp.(°C)	-60~-86	Shelves	3
Ambient Temp. (°C)	10-43	Noise(dBA)	≤53
Alarm System	High/low temperature, sensor failure, door ajar, power failure(72h), low battery		
Optional Accessories	Remote alarm system, cold storage racks & boxes, chart recorder		

## ST-MDF-86V838(VIP)



Model	ST-MDF-86V838	Capacity(L)	688
External Size(mm)	1160 X 1172 X 1907	N.W./G.W.(kg)	314.5/391
Inner Size(mm)	894 X 718 X 1310		
Package Size(mm)	1212 X 1098 X 2185	Insulation	PURF
Voltage	220V; 50/60Hz	External Material	Color sprayed steel
Power(W)	1540	Inner Material	304 Stainless steel
Setting Temp.(°C)	-60~-86	Shelves	3
Ambient Temp. (°C)	10-43	Noise(dBA)	≤53
Alarm System	High/low temperature, sensor failure, door ajar, power failure(72h), low battery		
Optional Accessories	Remote alarm system, cold storage racks & boxes, chart recorder		

## ST-MDF-86V936



Model	ST-MDF-86V936	Capacity(L)	936
External Size(mm)	1240 X 1000 X 2000	N.W./G.W.(kg)	340.5/416.5
Inner Size(mm)	1000 X 700 X 1335		
Package Size(mm)	1350 X 1135 X 2140	Insulation	PURF
Voltage	220V; 50/60Hz	External Material	Color sprayed steel
Power(W)	1600	Inner Material	304 Stainless steel
Setting Temp.(°C)	-40~-86	Shelves	3
Ambient Temp. (°C)	10-32	Noise(dBA)	≤53
Alarm System	High/low temperature, sensor failure, door ajar, power failure(72h), low battery		
Optional Accessories	Remote alarm system, cold storage racks & boxes, chart recorder		

## ST-MDF-86H105



Model	ST-MDF-86H105	Capacity(L)	105
External Size(mm)	680 X 580 X 935	N.W./G.W.(kg)	62/69.5
Inner Size(mm)	480 X 380 X 665		
Package Size(mm)	770 X 705 X 1090	Insulation	PURF
Voltage	220V; 50/60Hz	External Material	Color sprayed steel
Power(W)	280	Inner Material	304 Stainless steel
Setting Temp.(°C)	-40~-86	Compressor/QTY	Brand/1
Ambient Temp. (°C)	10-32	Noise(dBA)	≤52
Alarm System	High/low temperature, sensor		
Optional Accessories	Remote alarm system, cold storage racks & boxes		

## ST-MDF-86H485



Model	ST-MDF-86H485	Capacity(L)	485
External Size(mm)	1695 X 805 X 955	N.W./G.W.(kg)	124/138.5
Inner Size(mm)	1455 X 525 X 665		
Package Size(mm)	1775 X 900 X 1050	Insulation	PURF
Voltage	220V; 50/60Hz	External Material	Color sprayed steel
Power(W)	485	Inner Material	304 Stainless steel
Setting Temp.(°C)	-40~-86	Compressor/QTY	Brand/1
Ambient Temp. (°C)	10-32	Noise(dBA)	≤54.4
Alarm System	High/low temperature, sensor		
Optional Accessories	Remote alarm system, cold storage racks & boxes		

## ST-MDF-60V50



Model	ST-MDF-60V50	Capacity(L)	50
External Size(mm)	540 X 620 X 1115	N.W./G.W.(kg)	75/95
Inner Size(mm)	361 X 429 X 406		
Package Size(mm)	690 X 760 X 1170	Insulation	PURF
Voltage	220V; 50/60Hz	External Material	Color sprayed steel
Power(W)	420	Inner Material	304 Stainless steel
Setting Temp.(°C)	-25~-60	Shelves	1
Ambient Temp. (°C)	10-32	Noise(dBA)	≤49.5
Alarm System	High/low temperature, sensor		
Optional Accessories	Remote alarm system, cold storage racks & boxes, chart recorder		

## ST-MDF-60V340



Model	ST-MDF-60V340	Capacity(L)	340
External Size(mm)	854 X 1072 X 1907	N.W./G.W.(kg)	249.5/302
Inner Size(mm)	488 X 607 X 1140		
Package Size(mm)	935 X 1070 X 2135	Insulation	PURF
Voltage	220V; 50/60Hz	External Material	Color sprayed steel
Power(W)	1200	Inner Material	304 Stainless steel
Setting Temp.(°C)	-25~-60	Shelves	3
Ambient Temp. (°C)	10-43	Noise(dBA)	≤50
Alarm System	High/low temperature, sensor		
Optional Accessories	Remote alarm system, cold storage racks & boxes, chart recorder		

## ST-MDF-60H105



Model	ST-MDF-60H105	Capacity(L)	105
External Size(mm)	680 X 580 X 935	N.W./G.W.(kg)	62/69.5
Inner Size(mm)	480 X 380 X 665		
Package Size(mm)	770 X 705 X 1090	Insulation	PURF
Voltage	220V; 50/60Hz	External Material	Color sprayed steel
Power(W)	280	Inner Material	304 Stainless steel
Setting Temp.(°C)	-25~-60	Compressor/QTY	Brand/1
Ambient Temp. (°C)	10-32	Noise(dBA)	≤53
Alarm System	High/low temperature, sensor		
Optional Accessories	Remote alarm system, cold storage racks & boxes		

## ST-MDF-60H485



Model	ST-MDF-60H485	Capacity(L)	485
External Size(mm)	1695 X 805 X 955	N.W./G.W.(kg)	124/138.5
Inner Size(mm)	1455 X 525 X 665		
Package Size(mm)	1775 X 900 X 1050	Insulation	PURF
Voltage	220V; 50/60Hz	External Material	Color sprayed steel
Power(W)	485	Inner Material	304 Stainless steel
Setting Temp.(°C)	-25~-60	Compressor/QTY	Brand/1
Ambient Temp. (°C)	10-32	Noise(dBA)	≤53
Alarm System	High/low temperature, sensor		
Optional Accessories	Remote alarm system, cold storage racks & boxes		